STORZ

KARL STORZ—ENDOSKOPE

201123 20 Kaltlicht-Fontäne HALOGEN 150

201123 20 Cold light fountain HALOGEN 150

201123 20 Fuente de luz fría HALOGEN 150

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Instruction Manual

Manual de instrucciones

	20 1123 20	Kaltlicht-Fontäne HALOGEN 150
	20 1123 20	Cold light fountain HALOGEN 150
	20 1123 20	Fuente de luz fría HALOGEN 150

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**Important information
for users of KARL STORZ
instruments**

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Thank you for your expression of confidence in the KARL STORZ brand name. Like all of our other products, this product is the result of years of experience and great care in manufacture. You and your organization have decided in favor of a modern, high-quality item of equipment from KARL STORZ.

This instruction manual is intended to serve as an aid in the proper setup, installation, and operation of the HALOGEN 150 cold light fountain. All essential details of the equipment and all actions required on your part are clearly presented and explained. We thus ask that you read this manual carefully before proceeding to work with the equipment. Insert this manual in its protective wallet and keep it available for ready reference in a convenient and conspicuous location near the equipment.

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
© All product pictures, photos and descriptions and texts are the intellectual property of KARL STORZ GmbH & Co. KG
Utilization and duplication by third parties have to be expressly authorized by KARL STORZ GmbH & Co. KG
All rights reserved.

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Caution: Federal (USA) law restricts this device to sale by or on the order of a physician.

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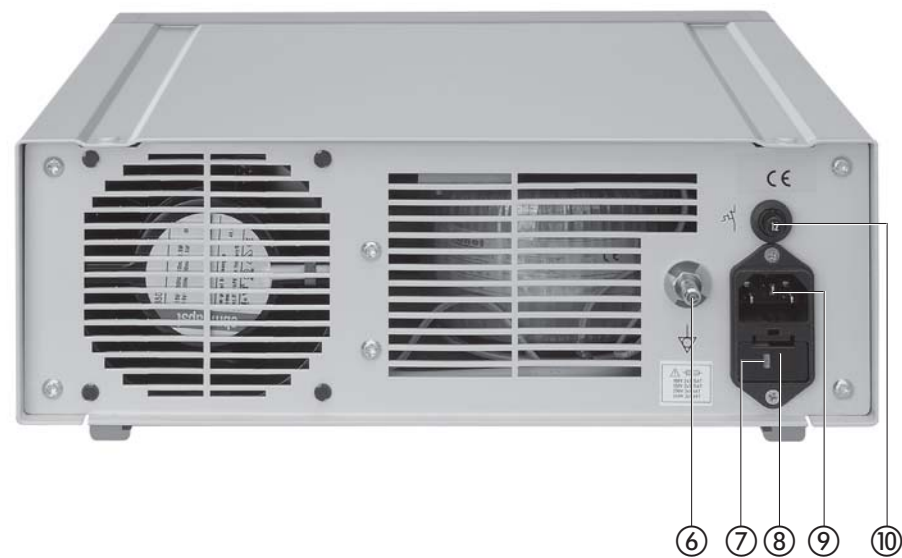**Controls, displays,
connectors and
their uses**

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
 Read the instructions carefully before
operating the equipment!

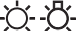 Luminous intensity switch

 Light exit

 Applied part of type CF

 Potential equalization connector

 Lamp fuse

 **Warning!** Burn hazard!

 Electronic information product
pollution control (China RoHS)

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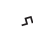
 **DANGER:** Risk of explosion if used
in the presence of flammable
anesthetics.

 **CAUTION:** To reduce the risk of
electrical shock, do not remove
cover. Refer servicing to qualified
service personnel.

Keep out of reach of patients.
Do not store liquids on or above the
unit.

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Photographs of the unit	A 1
Controls, displays, connectors, and their uses	A 2
Symbols employed	A 3

General information

Device description	2
--------------------------	---

Safety instructions

Warnings and cautions	3
Intended use	6
User qualification	7
Safety precautions at the site of installation	7
Safety precautions when operating the unit	7

Installation and operating instructions

Unpacking the equipment	8
Basic equipment	8
Installation	8
Operating the cold light fountain	9
Adaptor for light cables of other manufacturers	12

Maintenance

Lamp replacement	13
Fuse replacement	15
Cleaning and care of the unit	16
Maintenance	16
Servicing and repair	16
Repair program	17
Limitation of liability	18
Manufacturer's warranty	18

Technical description

Troubleshooting	19
Technical data	21
Technical documentation	22
General circuit diagram	23

Spare parts, recommended accessories

List of spare parts	24
Accessories	24

Appendix

Cleaning agents and disinfectants	30
Electromagnetic Compatibility (EMC) Information	32

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Device description

The cold light fountain HALOGEN 150 is suitable for all types of endoscopic examination. Color temperature is approx. 3.400 K.

A high-quality condenser system ensures high luminous efficiency while modern design guarantees easy operation and unproblematic servicing.

Adaptors for light cables of other manufacturers are available as accessories.

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Warnings and cautions

Please read this manual and follow its instructions carefully. The words **Warning**, **Caution**, and **Note** convey special meanings. Wherever they are used in this manual, they should be carefully reviewed to ensure the safe and effective operation of this product. To make the words, **Warning** and **Caution** stand out more clearly, they are accompanied by a pictogram.

Warning: A Warning indicates that the personal safety of the patient or physician may be involved. Disregarding a Warning could result in injury to the patient or physician.

Caution: A Caution indicates that particular service procedures or precautions must be followed to avoid possible damage to the product.

Note: A Note indicates special information about operating the product, or clarifies important information.

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Warning: Read this instruction manual thoroughly and be familiar with its contents prior to using this equipment. Before using the unit, read the following safety instructions carefully to avoid putting your patients, personnel or yourself at risk.

Note: Any damage to the unit resulting from incorrect operation is not covered by the manufacturer's warranty.

Warning: Test this equipment prior to each surgical procedure to ensure that it functions correctly.

Warning: Risk of explosion if used in the presence of flammable anesthetics.

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Safety instructions

Warnings and cautions

Warning: The electrical installations in the operating room in which the unit is connected and operated must comply with the applicable IEC standards.

Warning: Keep out of reach of patients.

Warning: The instructions and interface specifications for medical devices and/or system components used in combination must be observed precisely.

Warning: Grounding reliability can only be achieved when the equipment is connected to a properly installed "Hospital Only" or "Hospital Grade" receptacle (i.e., approved for use in an operating room environment). Routinely inspect electrical plug and cord. Do not use if inspection reveals damage.

Warning: Always unplug the unit before doing any maintenance work on it (e. g. cleaning).

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Warning: Danger of electric shock. Do not open the unit. Refer servicing only to qualified personnel. Any opening of the unit by unauthorized persons voids the guarantee. The only exception to this are activities described in this instruction manual.

Caution: Always adjust the cold light fountain to the minimum illumination intensity necessary to achieve optimum illumination of the endoscopic scene, either by direct vision or when coupled to a video camera.

Caution: Avoid allowing fluids to enter the unit. Do not store liquids on or above the unit.

Caution: Before connecting the unit to the electrical supply lines, verify that the line voltage stated on the unit's identification plate agrees with that of the electrical supply lines to be used.

Caution: Use only fuses of the correct rating.

Caution: Do not touch glass bulb with your bare fingers, use the protective cover. Touching the lamp with bare fingers may result in starting difficulty.

Caution: Federal law restricts this device to sale by or on the order of a physician (USA only).

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Warning: High energy radiated light through endoscopes may give rise to high temperatures in front of the light outlet and at the tip of the endoscope.

Warning: Never leave an illuminated light cable with or without an endoscope connected, on or near a patient; the risk of patient burns can result.

Caution: Do not place connected light cables on combustible materials such as textiles (drapes) or in the direct vicinity of swabs, etc. soaked with combustible fluids, such as disinfectant. These could catch fire from heat development at the open end of the light cable.

Therefore always attach the open end of a connected light cable to the provided retainer ④.

Note: Lamps do not have an indefinite service life, so it is possible that one may fail during use. For this reason a spare lamp system or cold light source should always be kept available during therapeutic endoscopic interventions.

Note: This unit has been marked in accordance with the European Directive 2002/96/EC on waste electrical and electronic equipment (WEEE). At the end of its useful operating life, dispose of the unit as electronic scrap.

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Before using the unit on the patient it is imperative that you be acquainted with how the unit operates and is controlled.

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Intended use

The cold light fountain HALOGEN 150 Model 20 1123 20 is used as a light source in conjunction with rigid and flexible endoscopes and other instruments that are designed for use with a cold light fountain.

Use of the unit in fields other than those indicated above is not allowed for safety reasons.

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The cold light fountain HALOGEN 150 may only be used with accessories which have been designated by KARL STORZ as suitable for the unit or the safe use of which is proven.

Note: KARL STORZ light cables are recommended for use with this unit. Some manufacturer's cables may be damaged when this unit is operated at higher settings.

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Unauthorized conversions or modifications to the unit are not allowed for safety reasons.

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
User qualification

The cold light fountain HALOGEN 150 may only be used by physicians and medical assistants who have a corresponding specialized qualification and who have been instructed in use of the unit.

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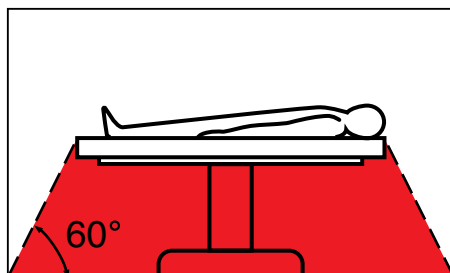Safety precautions at the site of installation

The unit may only be used in medical rooms whose electrical systems have been installed in accordance with applicable national regulations. It is not intended for use in hazardous zones.

This means, for example, that when using easily combustible and explosive inhalation anesthetics or mixtures thereof, the unit must not be operated inside the hazard zone shown in the diagram. This also applies for easily combustible and explosive chemicals, e.g. skin disinfectants and fast-acting surface disinfectants.

The unit is equipped with a connector for attaching a ground line. It should be connected before power is applied to the unit in compliance with current national regulations.

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Safety precautions when operating the unit

It is the user's responsibility to make sure the unit is safe and operates properly before using the unit.

During treatment using the cold light fountain HALOGEN 150 the patient must be treated and kept under observation with the usual medical care. This includes keeping a check on the progress of treatment, as well as monitoring sterile application conditions where required by the type of intervention.

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Unpacking the equipment

Carefully unpack the cold light fountain and remove the unit and its accessories from their packaging.

Check for missing items and evidence of shipping damage.

File any complaints with the manufacturer or supplier immediately. If possible, retain the original packing materials for later use; these can come in handy if the unit has to be transported.

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
Basic equipment

- 1 Cold light fountain HALOGEN 150
- 1 Power cord 400 A
- 1 Instruction manual
- 1 Wallet for instruction manual

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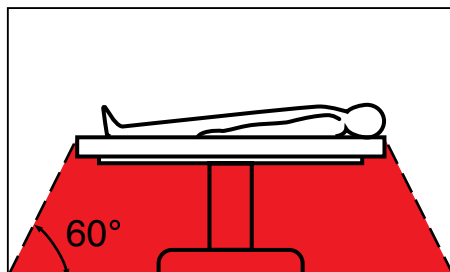Installation

Note: The cold light fountain HALOGEN 150 may be used only in medical facilities having electrical installations conforming to the applicable national, state, and local electrical codes.

Warning: This unit is not intended for use in hazardous zones. Do not operate the unit within the hazard zone shown in the diagram while explosive anesthetic gases are in use.

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Set the unit on a flat surface.

Caution: Do not cover the air vent. Danger of overheating. Make certain there is free air circulation towards the top and to the side (distance of rear panel to wall: 10 cm). If the unit overheats it will shut down automatically. Turn off the main power switch ① and let the light source cool down (see error list).

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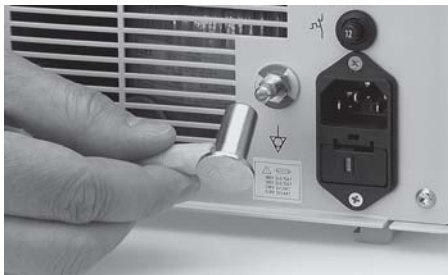
The unit is equipped with a connector ⑥ for attaching a ground line. The unit's ground line should be installed by a qualified electrician.

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
Operating the cold light fountain Setting of the system voltage

Remove the power fuse holder ⑧ (see section “Fuse replacement”, page 15) and take out the voltage selection module.

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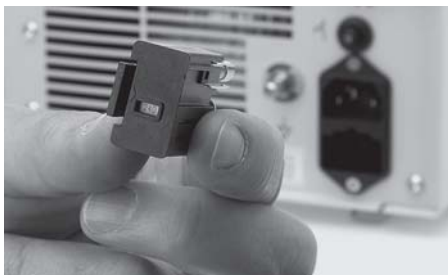Insert the voltage selection module into the power fuse holder in such a way that the system voltage displayed in the viewing window on the back of the power fuse holder corresponds to the local power supply.

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Check the inserted power fuses using the following table:

100 VAC	2 x T 3.15 AL250V
120 VAC	2 x T 3.15 AL250V
230 VAC	2 x T 1.6 AL250V
240 VAC	2 x T 1.6 AL250V

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Connect power cord. Insert power cord into power cord receptacle ⑨ as far as it will go.

Warning: The unit should only be connected to or disconnected from the power supply outside potentially explosive locations.

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Plug the light cable into socket ⑤ as far as it will go. To disconnect the light cable pull on the connector, never on the cable itself.

Warning: Danger of glare. Never look into the open end of a light cable.

Warning: Danger of burns. Contact with the open end of a light cable connected to a cold light fountain can possibly lead to burns.

Warning: Light cable socket may be very hot; particularly if the temperature switch has been triggered and the unit has switched off.

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Connect the light cable to the endoscope by a quarter turn of the knurled collar on the screw base.

Note: KARL STORZ light cables are recommended for use with this unit. Some manufacturer's cables may be damaged when this unit is operated at higher settings.

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If you do not plan to conduct the light cable to the endoscope right away, plug the loose end into the dummy socket ④.

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Note: The entrance and exit faces of the prepared light cable and endoscope must be kept clean if optimal light delivery is to be maintained. Remove heavy soiling with a lint-free disposable cloth, disinfect the entrance and exit faces with a cotton tip applicator soaked in alcohol and dry.

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Switch on the unit. The power switch ① and the operating pilot lamp ③ are lit.

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Note: The light exit ⑤ is provided with an anti-glare flap which does not allow any direct outlet of light. Light only appears when the light cable is attached.

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Adjust the required luminous intensity with the intensity switch ② (3 levels).

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Note: Frequently switching the cold light fountain on and off causes the halogen lamp to wear out more quickly. For relatively brief interruptions, set the luminous intensity switch ② to Level 1.

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Adaptor for light cables of other manufacturers

If necessary light exit ⑤ can be unscrewed and replaced by adaptors for light cables of other manufacturers (see section “Spare parts, recommended accessories”, page 29).

To do so, unscrew light exit ⑤ by rotating counterclockwise.

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Place the adaptor on the screw thread and tighten up by rotating clockwise.

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Lamp replacement

Switch off the unit. Disconnect the power plug.

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Open the housing cover. For this purpose, first undo the 2 thumbscrews at the bottom of the unit.

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Turn the unit back up the right way and slide the housing cover towards the back.

Note: Slide back the housing cover just far enough that the lamp compartment is easily accessible.

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
Remove the lid from the lamp compartment. To do this, loosen the thumbscrew.

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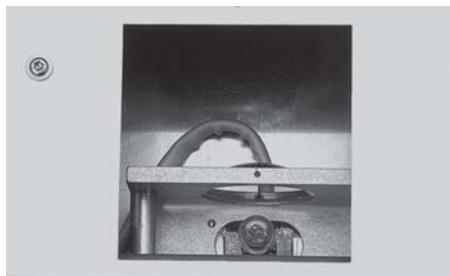Pull out the defective halogen lamp.

Warning: Danger of burns. Lamp may be hot.


Caution: Do not touch the glass of the new bulb with your fingers.

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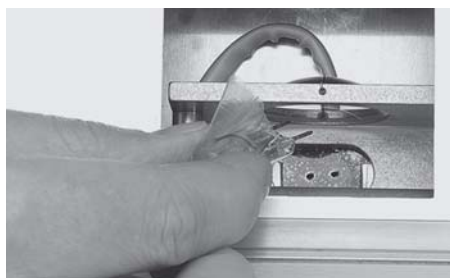Note: Do not touch glass bulb with your bare fingers. The glass bulb should be free of grease and other residue build-up. (If necessary, wipe the bulb with a cotton swab moistened with alcohol before operation).

Please use the protective cloth provided with the new lamp when inserting the lamp. Insert new halogen lamp (15 V/150 W) into the contact socket until firmly in place. The lamp is automatically adjusted.

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Reinstall the lamp compartment lid. Close the housing and tighten the thumbscrews. Reconnect the power cord. Carry out function test.

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Fuse replacement

Switch off the unit and remove the power cord.

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Remove the power fuse holder ⑧ with a suitable tool.

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! Caution: Only use fuses of the correct rating.

Insert new fuses of the appropriate rating.

	100/120 VAC	230/240 VAC
Power fuses	2 x T 3.15 AL	2 x T 1.6 AL

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Refit the power fuse holder ⑧.
Reconnect the power cord.
Test the instrument for proper operation.

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Cleaning and care of the unit

Warning: Always pull out power plug before cleaning! Make sure the lamp has cooled down.

Caution: Avoid allowing fluids to enter the unit.

The exterior of the unit can be surface disinfected with a disposable cloth soaked in disinfectant (preferably alcohol-free).

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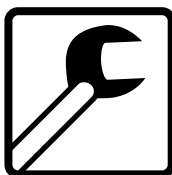
Maintenance

Preventive maintenance is not essential. Regular maintenance can, however, contribute to identifying potential problems before they become serious, thus enhancing the instrument's reliability and extending its useful operating life. Maintenance services can be obtained from your local representative or from the manufacturer.

Regardless of the accident prevention regulations or testing intervals for medical instruments prescribed in different countries, we recommend a functional or safety test of the unit at least once a year. For detailed instructions please refer to the version of the Service Manual currently valid, catalog number SV 3515 (English version).

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Servicing and repair

Defective items of equipment must be serviced and repaired exclusively by persons authorized by KARL STORZ; all repair work must employ original KARL STORZ parts only.

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Disposal

This unit has been marked in accordance with the European Directive 2002/96/EC on waste electrical and electronic equipment (WEEE). At the end of its useful operating life, dispose of the unit as electronic scrap.

Please ask either KARL STORZ GmbH & Co. KG, a KARL STORZ subsidiary or your specialist dealer for information on your local collection point. Within the scope of application of this Directive, KARL STORZ GmbH & Co. KG is responsible for the proper disposal of this unit.

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Repair program

Usually to bridge the repair period, you will receive a unit on loan which you then return to KARL STORZ as soon as you receive the repaired unit.

In Germany you can refer repairs direct to

KARL STORZ GmbH & Co. KG
Abt. Reparaturservice
Dr. Karl-Storz-Str. 34
78532 Tuttlingen.

In other countries please contact your local KARL STORZ branch or authorized dealer.

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Important information

In order to prevent infection, it is strictly forbidden to ship contaminated medical devices. All medical devices must be decontaminated on site to avoid contact and aerogenous infections (among staff). We reserve the right to return contaminated instruments/equipment to the sender.

Repairs, changes or expansions which are not carried out by KARL STORZ or by experts authorized by KARL STORZ will invalidate all guarantee rights.

KARL STORZ gives no guarantee on the correct functioning of equipment or instruments which have been repaired by unauthorized third parties.

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Limitation of liability

KARL STORZ shall be liable for failure or deterioration in the safe operation, operational reliability, and performance of this equipment only subject to the following conditions:

- all assembly operations, system expansions, readjustments, modifications, or repairs have been performed by a person or persons duly authorized by KARL STORZ;
- all electrical installations in the room where the unit is connected and operated comply with the applicable laws and standards; and
- the unit has been used in accordance with its operating instructions at all times.

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Manufacturer's warranty

For a period of **two** years after delivery to the end-user, we shall replace free of charge equipment that can be proven to have defects in material or workmanship. In so doing we cannot accept to bear the cost of transportation or the risk of shipment. The warranty referred to in our Standard Conditions of Business shall apply. Please fill out the attached warranty card on the last page and return it as soon as possible to:

KARL STORZ GmbH & Co. KG
Postfach 230
78503 Tuttlingen
Germany

Note: The unit must always be sent to the following address (also during the warranty period; where applicable, with warranty card):

KARL STORZ GmbH & Co. KG
Abt. Reparaturservice
Dr. Karl-Storz-Str. 34
78532 Tuttlingen
Germany

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Opening the equipment or performance of any repairs or modifications to the equipment by unauthorized persons shall relieve us of any liability for its performance. Any such opening, repair, or modification performed during the warranty period shall void all warranty.

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Troubleshooting

Warning: Always unplug the unit before all maintenance work.

Symptom:
Complete failure of the unit.

Possible causes:

- No power from the power line.
- Defective power fuse.

Remedy:

- Check that there is electricity to the wall outlet.
- Change fuses as described in the instruction manual. Make sure you use the correct type of fuse.

Symptom:
No light output, operating pilot lamp ③ is off. Power switch is lit (on).

Possible causes:

- Lamp defect.
- Temperature safety device triggered.

Remedy:

- Insert new lamp.
- Disconnect unit from power and allow to cool down. Switch on again after approx. 5 minutes. Make certain there is sufficient ventilation.

Symptom:
Too little light.

Possible causes:

- Luminous intensity switch is set too low (Stage 1 or 2).
- Ends of light cable and/or endoscope are soiled.
- Light cable is defective.
- Dimmer defective.

Remedy:

- Set the luminous intensity switch to a higher stage (2 or 3)
- Clean ends of light cables, and light entrances and exits on the endoscope.
- Replace light cable.
- Call service.

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Symptom:

No light output in spite of changing the bulb.

Possible causes:

- Bulb inserted is defective.
- Automatic cutout ⑩ has operated.
- Bulb base is defective.

Remedy:

- Fit a new bulb.
- Activate the automatic cutout ⑩ by pressing.
- Have the unit repaired.

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Standard compliance (for 201123 20)

According to IEC 60601-1:

- Type of protection against electric shocks: Protection Class I
- Degree of protection against electric shocks: Applied part of type CF 

According to IEC 60601-1-2:

Please read the Electromagnetic Compatibility Information in the appendix (p. 32-43).

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Directive compliance (for 201123 20)

According to Medical Device Directive (MDD):

This medical product belongs to Class I

This medical product bears the CE mark in accordance with MDD 93/42/EEC.

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Technical documentation

On request, the manufacturer will provide those circuit diagrams, itemized parts listings, descriptions, sets of adjustment instructions, and other items of available documentation to suitably qualified user personnel duly authorized by the manufacturer for their use in repairing those components of the instrument that have been designated by their respective manufacturers as repairable.

Supply of such technical documentation relating to the instrument shall not be construed as constituting manufacturer's authorization of user's personnel, regardless of their level of technical training, to open or repair the instrument.

Explicitly exempted herefrom are those maintenance and repair operations described in this manual.

We reserve the right to make engineering modifications in the interest of promoting technological progress and generating performance improvements without obligation on our part to submit prior notice thereof.

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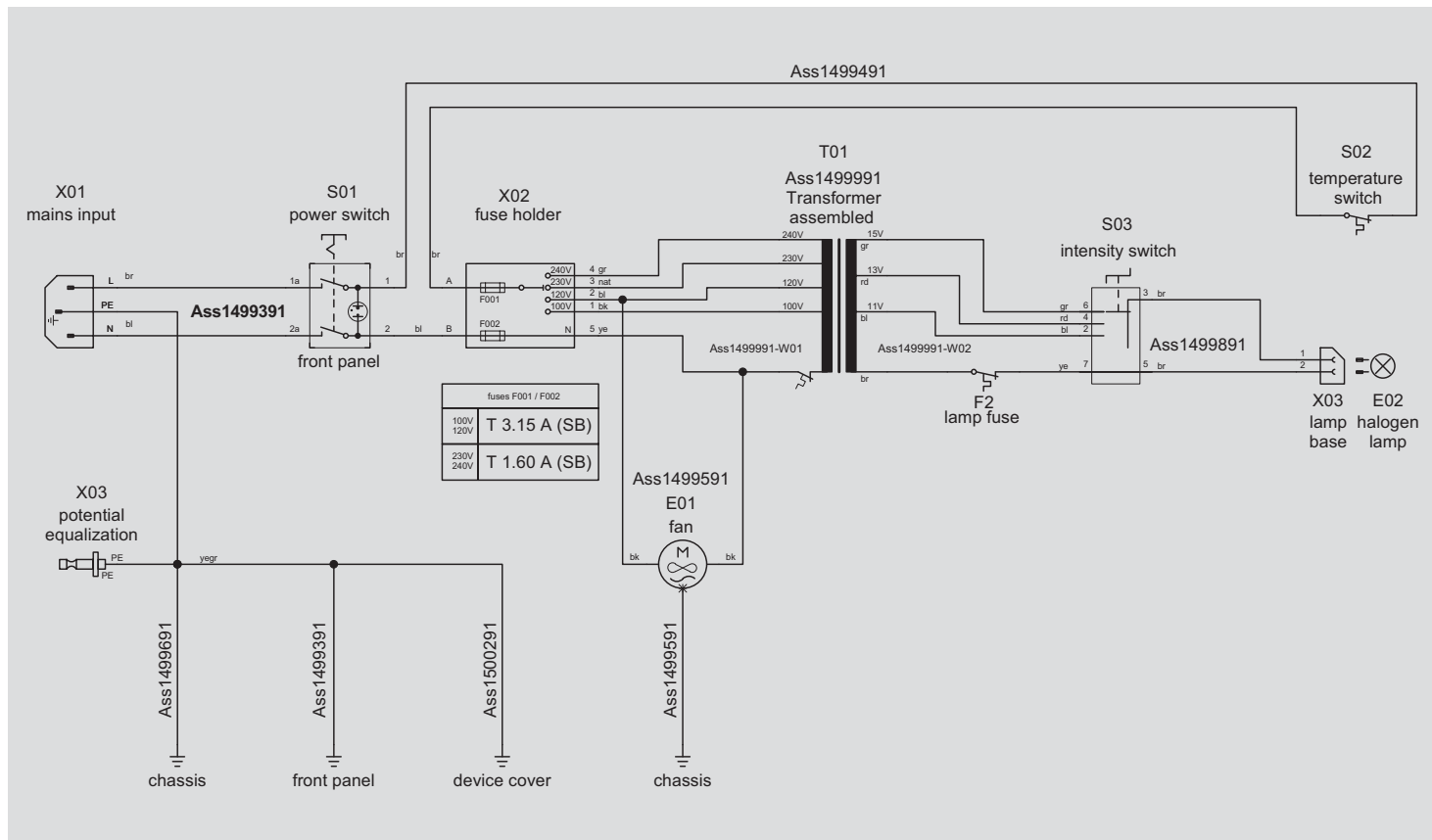
**Accessories
Fiberoptic light cable**

Fiberoptic light cables contain a large number of single glass fibers. The spectral transmission of a fiberoptic light cable is approximately the same for all of the wavelengths of visible light. It decreases, however, towards the blue end of the spectrum. This causes the light to take on a warmer shade – an effect which becomes more intensive the longer the light cable. Fiberoptic light cables are opaque to ultraviolet light.

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Cleaning, disinfection and sterilization of fiberoptic light cables

Note: Comprehensive information on cleaning light cables can be found in the 'Light cable' instructions.

1. **Manual (preliminary) cleaning:** The light cables may be cleaned using a cleaning agent/disinfectant which is specially approved for endoscopes. Wipe the outside of the light cables with a sponge or cloth soaked in cleaning solution.*
2. **Disinfection:** immerse in a disinfection solution (see Appendix, page 30, 31).*

* re 1./2. The light cable must always be rinsed until it is free of all residue after it has been cleaned or disinfected manually. Microbiologically pure/sterile water must be used for the final rinsing.

3. **Thermal decontamination by machine** (up to 93°C)

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4. **Manual final cleaning:** Clean the fiber surfaces and optical end faces with a soft cloth, sponge or cotton swab, soaked with 70% alcohol. Stubborn contaminants and deposits may be removed using cleaning paste (Art. no. 27661).
5. **Sterilization:** Steam sterilization at 134°C (+3°C) STERIS SYSTEM1® processor or gas or plasma sterilization.

Note: The instructions for 'Cleaning, Sterilization and Care of KARL STORZ Instruments' must be followed. There the procedures for cleaning, disinfection and sterilization are explained in detail.

Use a cleaning agent/disinfectant approved by KARL STORZ for cleaning/disinfection (see Appendix). You will also find the current list of approved products on our Web site (www.karlstorz.com).

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Fiberoptic light cable

Combination of light cable with endoscope

For safety reasons, the diameters of fiberoptic/ fluid* light cables and endoscopes should be used in the combination shown in the table opposite.
* see page 27

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Fluid light cables

In fluid light cables the light is not transmitted through glass fibers, but through a special liquid contained in the light cable. Fluid light cables are stiffer than glass fiber light cables and bending them too much will cause a severe loss of light. However, they do transmit a brighter light than fiberoptic light cables with the same diameter. For this reason, fluid light cables are ideally suited for endoscopic photography. Fluid light cables are permeable to ultraviolet light.


Caution: Damage to the elastomer coat will impair the light cable.

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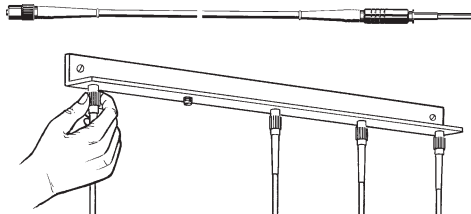Cleaning, disinfection and sterilization of fluid light cables

1. Manual (preliminary) cleaning: The fluid light cables may be cleaned using a cleaning agent/disinfectant which is specially approved for endoscopes. Remove heavy soiling with a moist sponge or lint-free disposable cloth.*

2. Disinfection: immerse in a disinfection solution (see Appendix, page 30, 31).*

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3. Machine decontamination

Caution: Carry out chemical disinfection by machine at **max. 65°C**.

Caution: When preparing light cables in a cleaning and disinfection machine, particular care is necessary in the choice of chemical agents and machine programs.

4. Manual final cleaning: Clean the fiber surfaces and optical end faces with a soft cloth, sponge or cotton swab, soaked with 70% alcohol.

Stubborn contaminants and deposits may be removed using cleaning paste (Art. no. 27661).

5. Sterilization: Gas sterilization or STERIS SYSTEM1® processor.

Caution: Fluid light cables must not be steam or plasma sterilized under any circumstances.

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More accessories are listed in our TELEPRESENCE catalogue.

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Typ

- 1 = Cleaner
2 = Disinfectant
3 = Suitable for ultrasound baths
A = Not suitable for flexible endoscopes (according to manufacturers' instructions)
(A)= Only for the flexible endoscopes of the new generation (year of manufacture > 2005)

Note: You will also find the current list of approved products on our Web site (www.karlstorz.com).

Note: When agents other than those named are used, KARL STORZ assumes no liability for any damage which may occur.

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Extract from the list of approved chemicals for the preparation of telescopes and endoscopic instrument sets
Manual preparation

The following instruments and telescopes produced by KARL STORZ must not be fully immersed in a liquid: HAMOU® telescopes with ocular focusing, magnifiers with focusing, IMPERATOR drill handles, prismatic light deflectors. For motors and their handles the appropriate manual provides detailed information.

The following instruments produced by KARL STORZ must not be placed in an ultrasound cleaning bath: rigid telescopes, flexible telescopes, other instruments with optical components made of glass.

When preparing and using the solutions, follow the solution manufacturer's instructions, paying close attention to concentration and exposure times. Regarding the microbiological effectiveness please ask the manufacturer. The following cleaning and disinfectant solutions have been approved regarding material compatibility:

Ecolab GmbH & Co. OHG	Sekusept Plus Sekusept Pulver Classic Ultrasemit AF Ultrassept Aktiv Ultrasseptin Classic	1, 2, 3 1, 2, 3 1, 2, 3 1, 2, 3 1, 2, 3
Esteer Pharma GmbH	Afid Afid plus Sporcid FF Polsept Holifa	1, 2 2 1, 2 1, 2
Fresenius AG	Letahdeido Darodor 4000 Liquido Darodor 9000	2 1 2
Holifa Polska Sp. zo.o	Darodor Sinaldehyd 2000	2
Indeba Indústria E Comercio LTDA	Aldasan 2000 Almyrol Desoform Lysoformin 3000	1, 2 1, 2 1, 2, 3 1, 2
Jose Collado S.A.	MedDis MediZyme Edisonite Super Mucadont-IS Mucadont-Zymaktiv Mucocit-1 neu	2 1 1, 3 2, 3 1, 3 1, 2, 3
Lysoform Dr. Rosemann GmbH	ID 50 A 20	2 1, 3, A
Medichem International	Wavacide Descoton Extra Descoton forte Perfektan Endo Perfektan Neu	2 1, 2 2 2, 3 1, 2
Merz Hygiene GmbH	Perfektan TB Gigasept Gigasept AF forte Gigasept FF	1, 2, 3, A 2, 3 1, 2 1, 2
METASYS Medizintechnik GmbH (formerly LONZA)	Gigasept Instru AF Gigasept Med forte Gigazyme Lysetol V	1, 2, 3, A 1, 2 1, 3 2
Orochemie	EnzyCare 2	2
Promagent AB	neodisher LM 2 neodisher mediclean neodisher mediclean forte neodisher medizym neodisher SeptoClean neodisher SeptoPreClean Aidul Plus Matrix	1, 3 1 1, 3 1, 3 1, A 1, 2 1 1
Schumacher, Dr. GmbH		
Schülke & Mayr GmbH		
Steris		
Weigert, Dr. GmbH & Co.		
Whiteley Industries PTY.Ltd		

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*** Caution:** The use of CAUSTIC SODA can cause changes to the surface of aluminium parts (including coated parts), plastics and soldered connections, and may reduce the service life of the instrument set. For **flexible endoscopes**, only immerse the examination sheath, **not** the housing/handle part.

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Typ

- 1 = Cleaner
- 2 = Disinfectant
- 3 = Suitable for ultrasound baths
- 4 = Neutralizer
- 5 = Rinse aid
- A = Not suitable for flexible endoscopes (according to manufacturers' instructions)

Note: You will also find the current list of approved products on our Web site (www.karlstorz.com).

Note: When agents other than those named are used, KARL STORZ assumes no liability for any damage which may occur.

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Machine preparation

The following instruments and telescopes produced by KARL STORZ must not be prepared in a cleaning and disinfecting machine: HAMOU® telescopes with ocular focusing, magnifiers with focusing, IMPERATOR drill handles, standard, tungsten carbide and diamond-tipped drills, prismatic light deflectors, and magnetic instrument rack. For motors and their handles the appropriate manual provides detailed information.

Selection of a method for cleaning and disinfection must be discussed with the machine manufacturer and the manufacturer of the chemical reagents. Only special methods must be used which have been verified for this purpose. Regarding the microbiological effectiveness please ask the manufacturer. Thermal disinfection is preferable. The following agents for machine cleaning and disinfection have been approved regarding material compatibility:

Schülke & Mayr GmbH	Thermosept alka clean	1
	Thermosept DK	2
	Thermosept ED	2
	Thermosept ER	1
	Thermosept RKN-zym	1
Steris	Instru-Klenz	1, 3
	neodisher DuoClean	1, A
Weigert, Dr. GmbH & Co.	neodisher FA	1
	neodisher FA forte	1, 3
	neodisher Oxivario	1
	neodisher mediclean	1
	neodisher mediclean forte	1, 3
	neodisher mediklar	5
	neodisher medizym	1
	neodisher SeptoClean	1, 2, A
	neodisher Septo DA	2
	neodisher Septo DN	2

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Appendix

Electromagnetic Compatibility (EMC) Information

WARNING: Medical electrical equipment needs special precautions regarding Electromagnetic Compatibility (EMC). Observe the EMC instructions in this appendix during installation and commissioning.

The HALOGEN 150 model 201123 20 corresponds to EN/IEC 60601-1-2 :2001 [CISPR 11 Class B] and therefore meets the EMC requirements of the Medical Device Directive (MDD) 93/42/EEC.

These limits are designed to provide reasonable protection against the typical electromagnetic interference to be expected in a medical environment. The HALOGEN 150 model 201123 20 is a Group 1 unit (as per CISPR 11). Group 1 contains all the 'equipment and systems which generate or use RF energy only for their internal functioning'.

NOTE: The tables and guidelines that are included in this appendix provide information to the customer or user that is essential in determining the suitability of the equipment or system for the electromagnetic environment of use, and in managing the electromagnetic environment of use to permit the equipment or system to perform its intended use without disturbing other equipment and systems or non-medical electrical equipment. If this equipment does cause harmful interference with other devices, the user is encouraged to try to correct the interference by one or more of the following measures:

- reorient or relocate the receiving device
- increase the separation between the equipment
- Connect the equipment into an outlet on a circuit different from that to which the other device(s) is connected.

If you have any further questions, please contact your local representative or our service department.

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Appendix

Electromagnetic Compatibility (EMC) Information

WARNING: The HALOGEN 150 model 201123 20 should not be used adjacent to or stacked with other equipment. If adjacent or stacked use is necessary, the equipment or system should be observed to verify normal operation in the configuration in which it will be used.

WARNING: The use of portable and mobile RF equipment may have an impact on this and other pieces of medical equipment.

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Table 201		
Guidance and manufacturer's declaration – electromagnetic emissions		
HALOGEN 150 model 201123 20 is intended for use in the electromagnetic environment specified below. The customer or user of the HALOGEN 150 model 201123 20 should ensure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	HALOGEN 150 model 201123 20 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause interference in nearby electronic equipment. The HALOGEN 150 model 201123 20 is suitable for use in all establishments (e.g. hospitals and doctors' offices) including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Complies	

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Table 202			
Guidance and manufacturer's declaration – electromagnetic immunity			
<p>HALOGEN 150 model 201123 20 is intended for use in the electromagnetic environment specified below. The customer or user of the HALOGEN 150 model 201123 20 should ensure that it is used in such an environment.</p>			
Immunity test	EN/IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic Discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	Complies ± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	Complies ± 2 kV for power supply lines ± 1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV differential mode ± 2 kV common mode	Complies ± 1 kV differential mode ± 2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<p><5% U_T * (>95% dip in U_T) for 0.5 cycle</p> <p>40% U_T (60% dip in U_T) for 5 cycles</p> <p>70% U_T (30% dip in U_T) for 25 cycles</p> <p><5% U_T (>95% dip in U_T) for 5 sec</p>	<p>Complies <5% U_T * (>95% dip in U_T) for 0.5 cycle</p> <p>40% U_T (60% dip in U_T) for 0.5 cycle</p> <p>See guidelines</p> <p>70% U_T (30% dip in U_T) for 1 cycle</p> <p>See guidelines</p> <p>Complies <5% U_T (>95% dip in U_T) for 5 sec</p>	<p>Mains power quality should be that of a typical commercial or hospital environment.</p> <p>If the user of the equipment or system requires continued operation during power cuts, short term power interruptions or if power supply is erratic, it is recommended that the equipment or system be powered from an uninterruptible power supply.</p>
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	Complies 3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
* Note: U_T is the a.c. mains voltage prior to application of the test level.			

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Table 204 Guidance and manufacturer's declaration – electromagnetic immunity – for equipment and systems that are not life-supporting			
HALOGEN 150 model 201123 20 is intended for use in the electromagnetic environment specified below. The customer or user of the HALOGEN 150 model 201123 20 should ensure that it is used in such an environment.			
Immunity test	EN/IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6 Radiated RF IEC 61000-4-3	$3 V_{\text{rms}}$ 150 kHz to < 80 MHz 3 V/m 80 MHz to 2.5 GHz	$3 V_{\text{rms}}$ 3 V/m	<p>Portable and mobile RF communications equipment should be used no closer to any part of the HALOGEN 150 model 201123 20, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance:</p> $d = [3.5/3]\sqrt{P} \quad 150 \text{ kHz to } < 80 \text{ MHz}$ $d = [3.5/3]\sqrt{P} \quad 80 \text{ MHz to } < 800 \text{ MHz}$ $d = [7/3]\sqrt{P} \quad 800 \text{ MHz to } 2.5 \text{ GHz}$ <p>Where P is the maximum output power rating of the transmitter in watts [W] according to the transmitter manufacturer and d is the recommended separation in meters [m].</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b.</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol:</p> 
Note: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection by structures, objects and people.			
<p>a Field strengths from fixed transmitters, such as base stations for radio [cellular/cordless] telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the unit is used exceeds the applicable RF compliance level above, the unit should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the unit.</p> <p>b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.</p>			

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Table 206

**Recommended separation distances between portable and mobile RF communications equipment and the
HALOGEN 150 model 201123 20**

The HALOGEN 150 model 201123 20 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or user of the unit can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the unit as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter [W]	Separation distance d [m] according to frequency of transmitter		
	150 kHz to < 80 MHz $d = \left[\frac{3.5}{3} \right] \sqrt{P}$	80 MHz to < 800 MHz $d = \left[\frac{3.5}{3} \right] \sqrt{P}$	800 MHz to 2.5 GHz $d = \left[\frac{7}{3} \right] \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance **d** in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where **P** is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

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WARRANTY

For a period of **two** years after delivery to the end-user, we shall replace free of charge equipment that can be proven to have defects in material or workmanship.

In doing so we cannot accept to bear the cost of transportation or the risk of shipment. The warranty referred to in our Standard Conditions of Business shall apply.

Complete warranty card upon purchase/delivery and return immediately to:

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Type of instrument: Serial no.:

Field of application:

To be filled out by instrument owner:
Return address/company stamp:

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KARL STORZ Endoskope
Regional Center for Endoscopy S.A.L.
St. Charles City Center, 5th Floor
Phoenicia Street, Mina Elhoshn
2020 0908 Beirut, Lebanon
Phone: +961 1 368181, Fax +961 1 365151

KARL STORZ Endoscopy South Africa, (Pty) Ltd.
8th Floor Convention Tower
Cnr Heerengracht & Coen Steytler, Foreshore
Cape Town 8001, South Africa
P.O.Box 6061, Roggebaai Cape Town,
8012, South Africa
Phone: +27 21 417 2600, Fax: +27 21 421 5103
E-Mail: info@karlstorz.co.za

TOO KARL STORZ ENDOSCOPY Kazakhstan
Khodjanova 17
050060 Almaty, Kazakhstan
Phone/Fax: +7 72 72 49 43 63, +7 72 72 49 41 00
e-mail: KarlStorzKazakhstan@itte.kz

KARL STORZ Endoscopy India Private Ltd.
C-126, Okhla Industrial Area, Phase-1
New Delhi 110020, India
Phone: +91 11 26815445-51,
Fax: +91 11 268129 86
E-Mail: karlstorz@vsnl.com

KARL STORZ GmbH & Co. KG
Resident Representative Office
80/33 (44/19) Dang Van Ngu,
F:10-Q, Phu Nhuan
Ho Chi Minh City, Vietnam
Phone: +848 991 8442, Fax: +848 844 0320

KARL STORZ Endoscopy (S) Pte. Ltd.
#05 – 08 San Centre, 171 Chin Swee Road
Singapore 0316, Singapore
Phone: +65 6376-1066, Fax: +65 6376-1068
E-Mail: karlstorz@pacific.net.sg

KARL STORZ Endoscopy Australia Pty. Ltd.
174 Parramatta Rd.
Camperdown, N.S.W. 2050, Australia
Phone: +61 02 8594 9150, Fax: +61 02 8594 9199
E-Mail: info@karlstorz.au

KARL STORZ Endoscopy China Ltd.
Hong Kong Representative Office
Unit 1601, Chinachem Exchange Square
1 Hoi Wan Street, Quarry Bay,
Hong Kong, People's Republic of China
Phone: +8 52 28 65 2411, Fax: +8 52 28 65 4114
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy China Ltd.
Beijing Representative Office
Room 610, China Life Tower
No. 6, Chaowai Street
Beijing, 100020, People's Republic of China
Phone: +86 10 8525 3725, Fax: +86 10 8525 3728
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy China Ltd.
Unit 3901-3904, Tower 1 Grand Gateway
No. 1 Hong Qiao Road
Shanghai, 200030, People's Republic of China
Phone: +86 21 6113-1188, Fax: +86 21 6113-1199
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy China Ltd.
Chengdu Representative Office
F-5, 24/F., Chuanxing Mansion,
No. 18 Renming Road South
Chengdu, Sichuan, 610016, People's Rep. of China
Phone: +86 28 8620-0175, Fax: +86 28 8620-0177
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy China Ltd.
Guangzhou Representative Office
Room 1119-20, Dongshan Plaza
69 Xianlie Road Middle, Dongshan District,
Guangzhou, Guangdong, 510095,
People's Rep. of China
Phone: +86 20 8732-1281, Fax: +86 20 8732-1286
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy Asia Marketing Pte Ltd
3791 Jalan Bukit Merah
06-11 e-Centre @ Redhill
Singapore 159471, Singapore
Tel. No. +65 63761066, Fax. No. +65 63761068
Email : infoasia@karlstorz.com.sg
Email : serviceasia@karlstorz.com.sg

KARL STORZ Endoscopy Singapore Sales Pte Ltd
3791 Jalan Bukit Merah
10-18 e-Centre @ Redhill
Singapore 159471
Tel. No. +65 65325548, Fax No. +65 65323832
Email : infosp@karlstorz.com.sg
Email : servicesgp@karlstorz.com.sg

KARL STORZ Endoscopy Japan K. K.
Bois Hongo Building 6Fl, 3-42-5 Hongo
Bunkyo-ku, Tokyo 113-0033, Japan
Phone: +81 3 5802-3966, Fax: +81 3 5802-3988
E-Mail: info@karlstorz.co.jp

www.karlstorz.com

KARL STORZ GmbH & Co. KG

Mittelstraße 8
78532 Tuttlingen

Postfach 230
78503 Tuttlingen
Germany

Telefon: +49 (0)7461 708-0
Telefax: +49 (0)7461 708-105
E-Mail: info@karlstorz.de
Web: www.karlstorz.com